
# Panasonic

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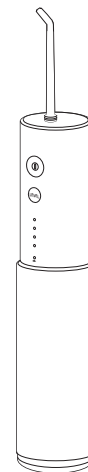## Operating Instructions (Household use) Oral Irrigator Model No. EW-DJ31

### Contents

<b>Safety precautions .....</b>	<b>4</b>	Frequently asked questions.....	15
Intended use.....	8	Troubleshooting.....	16
Parts identification .....	9	Battery life.....	18
Charging .....	9	Removing the built-in rechargeable	
How to use.....	11	battery.....	18
Points of care.....	12	Specifications .....	19
After use .....	13	Changing the nozzle (Sold separately)...	19

Thank you for purchasing this Panasonic product.

**Before operating this unit, please read these instructions completely and save them for future use.**

## Safety precautions

To reduce the risk of injury, loss of life, electric shock, fire, and damage to property, always observe the following safety precautions.

### Explanation of symbols

The following symbols are used to classify and describe the level of hazard, injury, and property damage caused when the denotation is disregarded and improper use is performed.

	<b>DANGER</b>	Denotes a potential hazard that will result in serious injury or death.
	<b>WARNING</b>	Denotes a potential hazard that could result in serious injury or death.
	<b>CAUTION</b>	Denotes a hazard that could result in minor injury or property damage.

The following symbols are used to classify and describe types of instructions to be observed.

	This symbol is used to alert users to a specific operating procedure that must not be performed.
	This symbol is used to alert users to a specific operating procedure that must be followed in order to operate the unit safely.

## **WARNING**

### ► Preventing accidents

**Do not store within the reach of children or infants. Do not let them use it.**

- Doing so may cause an accident or injury due to accidental ingestion of removable parts.

### ► Power supply

**Do not connect or disconnect the adaptor or the appliance plug to a household outlet with a wet hand.**

- Doing so may cause electric shock or injury.

**Do not immerse the AC adaptor in water or clean it with water.**

- Doing so may cause electric shock or fire due to a short circuit.

**Do not place the AC adaptor over or near water filled sink or bathtub.  
Do not use when the AC adaptor is damaged or when the fitting into the household outlet is loose.**

- **Do not damage, modify, forcefully bend, pull, twist, or bundle the cord. Also, do not place anything heavy on or pinch the cord.**  
- Doing so may cause electric shock or fire due to a short circuit.

## WARNING

**Do not use in a way exceeding the rating of the household outlet or the wiring.**

- Exceeding the rating by connecting too many power plugs to one household outlet may cause fire due to overheating.

**Do not use anything other than the supplied AC adaptor for any purpose. Also, do not use any other product with the AC adaptor.**

(See page 9.)

- Doing so may cause burn or fire due to a short circuit.

**Always ensure the appliance is operated on an electric power source matched to the rated voltage indicated on the AC adaptor. Fully insert the adaptor.**

- Failure to do so may cause fire or electric shock.

**The supply cord cannot be replaced. If the cord is damaged, the AC adaptor should be scrapped.**

- Failure to do so may cause an accident or injury.

### ► In case of an abnormality or malfunction

**Immediately stop using and remove the AC adaptor if there is an abnormality or malfunction.**

- Failure to do so may cause fire, electric shock or injury.

**<Abnormality or malfunction cases>**

- The main unit, adaptor or cord is deformed or abnormally hot.
- The main unit, adaptor or cord smells of burning.
- There is abnormal sound during use or charging of the main unit.
- Immediately request inspection or repair at an authorized service centre.

### ► This product

**This product has a built-in rechargeable battery. Do not throw into fire, apply heat, or charge, use, or leave in a high temperature environment.**

- Doing so may cause overheating, ignition, or explosion of the built-in rechargeable battery.

**Do not modify or repair.**

- Doing so may cause fire, electric shock, or injury.
- Contact an authorized service centre for repair (battery change, etc.).

**Never disassemble except when disposing of the product.**

- Doing so may cause fire, electric shock, or injury.

## WARNING

**This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.**

- Failure to do so may cause an accident or injury.

### ► Cleaning

**Always unplug the adaptor from a household outlet when cleaning the AC adaptor.**

- Failure to do so may cause electric shock or injury.

**Regularly clean the power plug and the appliance plug to prevent dust from accumulating.**

- Failure to do so may cause fire due to insulation failure caused by humidity. Disconnect the adaptor and wipe with a dry cloth.

## CAUTION

### ► Preventing damage to teeth and gums

**People who are unable to operate should not use this appliance.**

**People with weak oral senses should not use this appliance.**

- Doing so may cause injury or damage to teeth and gums.

**People who may have periodontal disease, have teeth being treated, or concerned about symptoms within their mouth should consult a dentist before use.**

- Failure to do so may cause injury or damage to teeth and gums.

### ► This product

**Never use for a purpose other than oral cleaning.**

- Doing so may cause an accident or injury when the stream of water is pointed toward the eyes, nose, ears, tongue, or throat.

**Do not fill the tank with water warmer than 40 °C.**

- Doing so may cause burns.

**Do not share nozzle with your family or other people.**

- Doing so may cause infection or inflammation.

## CAUTION

### ► Note the following precautions

**Do not allow metal objects or trash to adhere to the power plug or the appliance plug.**

- Doing so may cause electric shock or fire due to a short circuit.


**Do not drop or subject to shock.**

- Doing so may cause injury.

**Do not put your fingers on the nozzle storage unit or just hold by your hand the portion within the dotted line frame when storing.**

- Doing so may get your fingers/hand pinched.

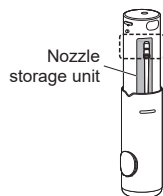**Do not wrap the cord around the adaptor when storing.**

- Doing so may cause the wire in the cord to break with the load, and may cause fire due to a short circuit.

**Disconnect the adaptor from the household outlet when not charging.**

- Failure to do so may cause electric shock or fire due to electrical leakage resulting from insulation deterioration.

**Disconnect the adaptor or the appliance plug by holding onto it instead of the cord.**

- Failure to do so may cause electric shock or fire due to a short circuit.

**After use, clean the main body, water tank and the nozzle and dry them. (See pages 13 to 15.)**

- Failure to do so may cause damage to your health due to propagation of mold and bacteria.

### Handling of the removed battery when disposing

## DANGER

**The rechargeable battery is exclusively for use with this appliance. Do not use the battery with other products.**

**Do not charge the battery after it has been removed from the product.**

**Do not throw into fire or apply heat.**

**Do not hit, disassemble, modify, or puncture with a nail.**

**Do not let the positive and negative terminals of the battery get in contact with each other through metallic objects.**

**Do not carry or store the battery together with metallic jewellery such as necklaces and hairpins.**

**Do not use or leave the battery where it will be exposed to high temperatures, such as under direct sunlight or near other sources of heat. Never peel off the tube.**

- Doing so may cause overheating, ignition, or explosion.

## WARNING

**After removing the rechargeable battery, do not keep it within reach of children and infants.**

- The battery will harm the body if accidentally swallowed.

If this happens, consult a physician immediately.

**If the battery fluid leaks out and comes in contact with the eyes, do not rub the eyes, and thoroughly rinse with clean water, such as tap water.**

- Failure to do so may cause eye injury.  
Consult a physician immediately.

## CAUTION

**If the battery fluid leaks out and comes in contact with the skin or clothes, rinse off with clean water, such as tap water.**

- Failure to do so may cause inflammation.

## Intended use

- Please use the appliance together with toothbrush. It is hard to remove tartar or dental calculus attached to the tooth surface just by a water jet, so a toothbrush is required to brush the teeth. We recommend rinsing the inside of your mouth with the water jet after brushing your teeth.
- When you first begin using the oral irrigator, you may experience slight bleeding of the gums even though


your gums may be healthy. This is because your gums are receiving stimulation by the oral irrigator for the first time. The bleeding should stop after one or two weeks of use. Bleeding that persists for more than two weeks may be indication of an abnormal condition of your gums, such as periodontal disease. In this case, stop using the oral irrigator and consult your dentist.

- Do not put salt, mouthwash or other chemicals into the water tank. Doing so may cause malfunction.
- Do not dry the main unit with a dryer or a fan heater. Doing so may cause malfunction or deformation of the parts.
- Store the appliance in a place with low humidity after use. Leaving the appliance in a bathroom may cause malfunction.
- Wrap the appliance in a towel or the like when carrying it with you. Otherwise, the water may seep out, thus wetting surrounding objects.
- Charge the appliance correctly according to these Operating Instructions. (See pages 9 and 10 "Charging".)
- This appliance contains battery that is only replaceable by skilled persons. Please contact an authorized service centre for the details of repair.
- The following symbol indicates that a specific detachable power supply unit is required for connecting the electrical appliance to the supply mains. The type reference of power supply unit is marked near the symbol.

## Parts identification
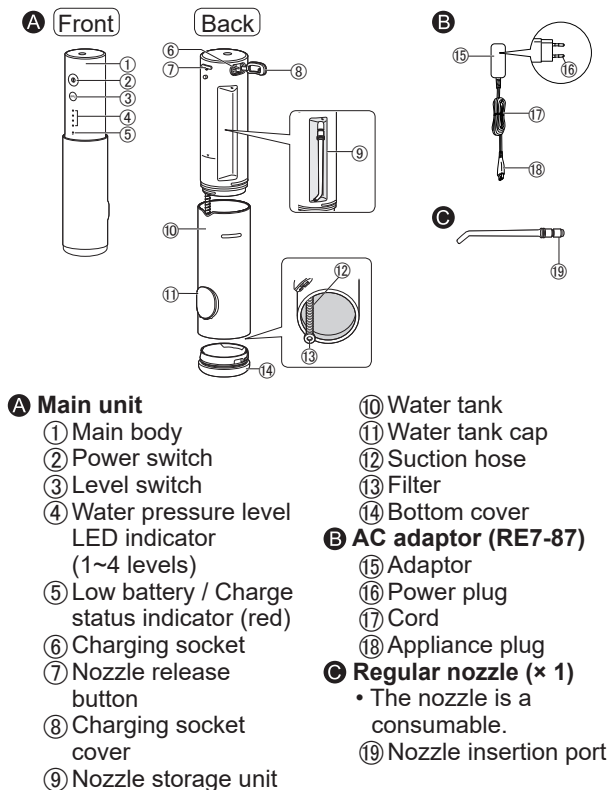
### Note


There may be moisture inside the main unit (water tank and suction hose) left from the distilled water used for product testing.

## Charging

- Wipe the charging socket dry before use if there are water drops in it.


### 1 Connect the AC adaptor.

- ① Open the charging socket cover, and then insert the appliance plug into the charging socket on the main body.
- ② Insert the adaptor into a household outlet.
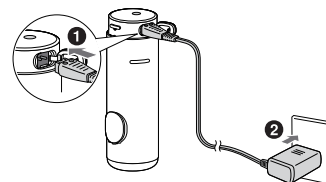
### ► About the charging reminder function

- During use, charging is required when the low battery / charge status indicator (red) flashes (once a second).
- Turning off the power switch "①" when the battery has little remaining charge, the low battery / charge status indicator will light off after 5 seconds.
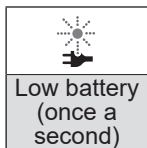
### ► Indicator display during charging


Charging in progress	Charging complete	Charging abnormal (Flashing in red twice a second)

(Lights on • Lights off ○) (Flashing ✨)

### ► How to check if charging is completed

- After charging is completed, if you disconnect the power plug and then plug in the power plug again to charge, (the low battery / charge status indicator) lights on and then lights off after about 5 seconds.

## 2 For safety and energy saving, disconnect the power plug after charging is completed.

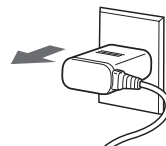- Disconnect the appliance plug from the charging socket, and then close the charging socket cover.
- To keep charging the battery will not affect battery performance.

### ► Operating time after a full charge approx. 10 minutes\*

(When the water pressure level is set to 4)

- \* It varies with the ambient temperature and usage conditions.
- Charging is completed after approx. 1 hour. (The charging time may vary with the remaining battery capacity.)


### Notes

- The charging time may differ when the appliance is used for the first time, or if it has not been used for more than 6 months. In addition, (the low battery / charge status indicator) may not light on in the first several minutes of charging, but it will light on as the charging proceeds.
- The recommended ambient temperature for charging is 10 °C - 35 °C. It may take longer to charge, or will flash rapidly and the battery may fail to charge properly in extremely high or low temperatures.

## How to use


### 1 Pull out the main body while rotating it, and then install the nozzle.

- ① Open the water tank cap.
  - If it is closed, the main body is hard to be pulled out.
- ② Pull the main body to the "▼" mark.
  - Be noted that the nozzle may be ejected from the nozzle storage unit.
- ③ Install the nozzle while holding the main body.
  - Install it into place until it clicks.

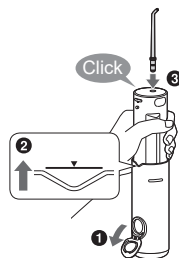
### 2 Fill the water tank with water.

- ① Add water while keeping the main unit level.
  - Do not fill beyond the full-water position.
- ② Close the water tank cap securely.
  - Always use new tap water or lukewarm water.
  - If you feel the water is cold, use lukewarm water with a temperature of 40 °C or less.

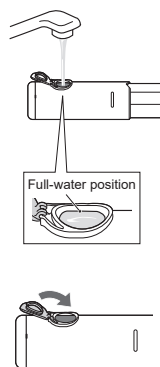
### ► If the water tank cap comes off

- Align the groove of the water tank cap with the protrusion of the water supply port, and press the water tank cap into place until it clicks.
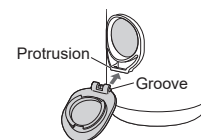
### 3 Press the level switch to select your desired water pressure level.

- After you press the level switch, the last used water pressure level is automatically selected. Then each time you press the level switch, the water pressure level is switched once in the following order: 1 (weak) → 2 → 3 → 4 (strong) → 1...

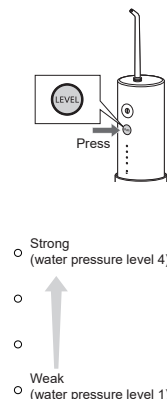
► **About the water pressure level memory function**

**The last selected water pressure level will be remembered, from which the appliance will start next time.**

- When the appliance has not been used for a long time or when the battery is used up, the water pressure level will return to 1.

## 4 **Clean inside your mouth with water jet.**

- ① Place the nozzle in your mouth.
  - ② Close your mouth gently.
  - ③ Turn on the power switch "①".
- You can switch the water pressure level by pressing the level switch, even while using the appliance.
  - If the main unit is tilted too much, the water does not come out from the nozzle.
- **A fully filled water tank provides approximately 35 seconds of use.**  
(When the water pressure level is set to 4)
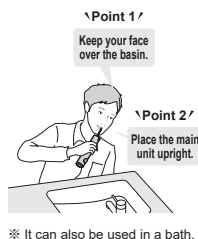
## Notes

- The appliance does not work when plugged into a household outlet.
- For the first use, turn on the power switch "①" to jet water so that you can check the water pressure level. (Please start using the appliance from water pressure level 1.)
- The appliance may not operate in an ambient temperature approximately 5 °C or lower.
- Do not push the tip of the nozzle too hard against the teeth or gums.
- If you feel that the irritation is strong, turn off the power switch "①" and use the appliance from water pressure level 1.
- Do not turn on the power switch "①" when the water tank is empty unless you are drying the moisture inside the main unit. (Doing so may cause malfunction.)

## Points of care

Change the way you use and move the unit depending on the purpose and the area you want to clean.

### Interdental cleaning

- Direct the water jet onto the space between the teeth.
- Direct water jet onto the front and rear of the teeth.

## Periodontal pocket cleaning

- Direct the water jet onto the space between the teeth and gums.
- Slowly shift the water jet along the gums.
- Clean the front and rear sides of the teeth as well as behind the molars.

### ► About tooth tartar

It is hard to remove tooth tartar or dental calculus attached to the tooth surface just by a water jet, so a toothbrush is required to brush the teeth. It is recommended to clean inside your mouth with the water jet after brushing your teeth.

### ► About bleeding gums

Initially bleeding may occur even with healthy gums. This is caused because the gums are not used to such stimulation. The bleeding will normally stop within 1 or 2 weeks of use. There may be a problem with the gums (periodontal disease, etc.) if the bleeding continues for over 2 weeks. In such cases, stop using the appliance and consult a dentist.

## After use


Please clean the appliance after each use.

- 1 To prevent the water from splashing about, turn off the power switch "①" while keeping the nozzle in your mouth.**


- 2 Empty the water tank and drain the remaining water in the main unit, for the sake of hygiene.**

- ①** Open the water tank cap, and pour out the remaining water.
  - Be sure to pour out the remaining water after each use.

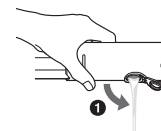
- ②** Turn on the power switch "①" to drain the remaining water in the main unit.
- ③** Turn off the power switch "①" when no more water comes out.
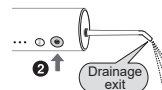
### 3 Remove the bottom cover, and wash the main body (as a whole) and inside of the tank.

- ① Remove the nozzle while pressing the nozzle release button.
- ② Remove the water tank from the main unit, and then take out the bottom cover.
  - **How to remove the bottom cover**
    - ① Rotate the bottom cover until the "I" mark on the water tank is aligned with the "I" mark on the bottom cover.
    - ② Remove the bottom cover.
- ③ Clean it with water.
  - When draining, do not hit the appliance on the washbasin or the like.

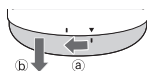
### 4 Store it after it is dry.

- ① Wipe away the water with a towel or the like and dry them in a well-ventilated shade.


### 2 Store the nozzle and the main unit.

- Attach the bottom cover onto the water tank.
- Do not get your fingers pinched in [ ] on the nozzle storage unit.
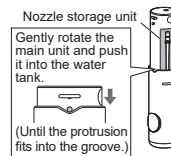
### ► How to attach the bottom cover

- ① Align the "I" mark on the water tank with the "I" mark on the bottom cover, and attach the bottom cover.
- ② Rotate the bottom cover until the "I" mark on the bottom cover is aligned with the "▼" mark on the water tank.
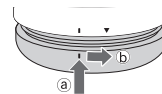
### ► If the rubber ring of the bottom cover comes off

- Keep the convex part of the rubber ring facing outward, and install the rubber ring onto the bottom cover.

## Notes


- Do not wash with lukewarm water exceeding 50 °C.
- Use neutral kitchen detergent and rinse well with water to ensure no detergent remains after cleaning.
- Never use thinners, benzine, or alcohol.  
(Doing so may cause malfunction, or cracking or discolouration of the parts.)
- The water tank and the bottom cover can be washed in a dish washer.

### ► When the water jet is poor

Remove the dirt from the filter (the tip of the suction hose).

#### Filter

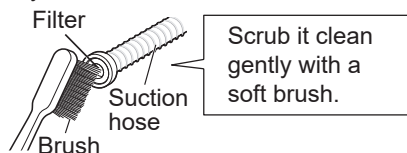Rub it gently with a soft brush.

## Frequently asked questions

Question		Answer
Battery	Will the battery be exhausted after a long period of disuse?	► Please charge it once every 6 months if you have not used it for over 6 months since the battery may deteriorate.
	Can the battery be charged before every use?	► This appliance uses a Lithium-ion battery. Therefore, the service life of the battery will not be affected even if it is charged after each use.

## Troubleshooting

	Problem	Possible cause	Action
While charging	The low battery / charge status indicator  does not light on.	<ul style="list-style-type: none"> <li>▶ The appliance was just purchased or it has not been used for a while.</li> <li>▶ The main unit is not properly plugged in the charging socket.</li> </ul>	<ul style="list-style-type: none"> <li>▶ Keep charging, and the indicator will light on after a while. (See pages 9 to 10.)</li> <li>▶ Check to see that the appliance plug is in full contact with the charging socket, and that the low battery / charge status indicator is on. (See pages 9 to 10.)</li> </ul>
	The low battery / charge status indicator flashes rapidly (twice a second).	<ul style="list-style-type: none"> <li>▶ The ambient temperature is too low or too high during charging.</li> </ul>	<ul style="list-style-type: none"> <li>▶ Please charge within the recommended temperature range of 10 °C to 35 °C.</li> <li>▶ If the indicator flashes even when you charge it within the recommended temperature range, contact your dealer or Service Centre for repair.</li> </ul>
	The appliance can be used for just several minutes even after a full charge.	<ul style="list-style-type: none"> <li>▶ The main unit is not properly plugged in the charging socket.</li> </ul>	<ul style="list-style-type: none"> <li>▶ Check to see that the appliance plug is in full contact with the charging socket, and that the low battery / charge status indicator is on. (See pages 9 to 10.)</li> </ul>
	The appliance does not work, but  does not flash as expected.	<ul style="list-style-type: none"> <li>▶ The battery has reached the end of its life (approx. 3 years).</li> </ul>	<ul style="list-style-type: none"> <li>▶ Have the battery replaced by your dealer or an authorized service centre.</li> </ul>

A diagram of a gas syringe. A label 'a' points to the stopper at the mouth of the syringe.

17

## Battery life

Approximately 3 years. The battery may have reached its operating life if the operating time decreases drastically, even after a full charge.

(The operating life of the battery may differ depending on usage or storage conditions.)

## Removing the built-in rechargeable battery


**Remove the built-in rechargeable battery before disposing of the appliance.**

The battery is to be disposed of safely.

Please make sure that the battery is disposed of at an officially designated location if there is one.

This figure must only be used when disposing of the appliance, and must not be used to repair it. If you dismantle the appliance yourself, it will no longer be waterproof, which may cause it to malfunction.

- Remove the appliance from the adaptor, press the power switch to turn on the power and then keep the power on until the battery is completely discharged.
- Use a screwdriver to perform steps ① to ⑦ to disassemble the main unit and remove the battery.
- Insulate the  $\oplus$  and  $\ominus$  terminals by applying tape to them. (To prevent short circuit of terminals).
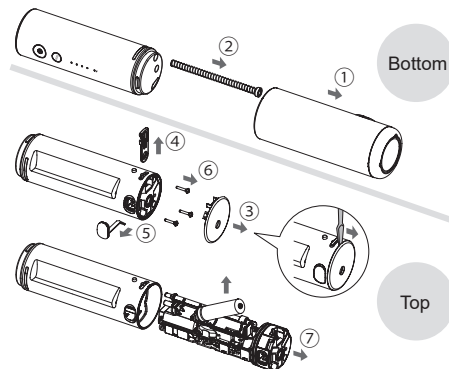
### For environmental protection and recycling of materials

This appliance contains a Lithium-ion battery.

Please make sure that the battery is disposed of at an officially assigned location, if there is one in your country.

## Specifications

Power source	See the name plate on the product (Automatic voltage conversion)
Charging time	Approx. 1 hour
Operating time	LEVEL 4: Approx. 10 minutes (at 20 °C when fully charged)
Tank volume	Approx. 150 mL (for approx. 35 seconds of operation when the water pressure level is set to 4)
Dimensions	Approx. 30.9 cm (or 23.4 cm when the nozzle is not installed) (H) x 5.2 cm (W) x 5.7 cm (D) In storage: Approx. 15.4 cm (H) x 5.2 cm (W) x 5.7 cm (D)
Mass	Main body: Approx. 240 g (Including nozzle)

This product is intended for household use only.

## Changing the nozzle (Sold separately)

### ► Regular nozzle

- If the tip of the nozzle gets deformed, please replace the nozzle with a new one.
- The nozzle should be replaced once every 6 months or so due to hygiene reasons, even if it is not deformed.

New nozzle (tip)

Deformed and discolored nozzle (tip)

Part number	Replacement Nozzle for EW-DJ31	
WEW0987	Regular nozzle	

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© Panasonic Corporation 2022
Panasonic Corporation