# Panasonic


## Operating Instructions

(Household use) Oral Irrigator

Model No. EW-DJ86

English.....3

Tiếng Việt ..... 19
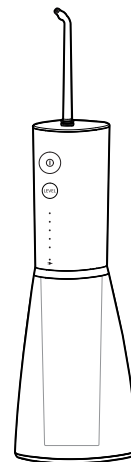
# Panasonic

## Operating Instructions

(Household use) Oral Irrigator

Model No. EW-DJ86

### Contents

|                                    |           |                                                            |           |
|------------------------------------|-----------|------------------------------------------------------------|-----------|
| <b>Safety precautions.....</b>     | <b>4</b>  | <b>Changing the nozzle<br/>(Sold separately).....</b>      | <b>13</b> |
| <b>Important information .....</b> | <b>6</b>  | <b>Troubleshooting .....</b>                               | <b>14</b> |
| <b>Parts identification .....</b>  | <b>7</b>  | <b>Battery life .....</b>                                  | <b>16</b> |
| <b>Charging.....</b>               | <b>8</b>  | <b>Removing the built-in rechargeable<br/>battery.....</b> | <b>16</b> |
| <b>How to use .....</b>            | <b>9</b>  | <b>Specifications .....</b>                                | <b>17</b> |
| <b>Points of care.....</b>         | <b>11</b> |                                                            |           |
| <b>Cleaning .....</b>              | <b>12</b> |                                                            |           |

Thank you for purchasing this Panasonic product.

Before operating this unit, please read these instructions completely and save them for future use.

## Safety precautions

To reduce the risk of injury, loss of life, electric shock, fire, and damage to property, always observe the following safety precautions.

### Explanation of symbols

The following symbols are used to classify and describe the level of hazard, injury, and property damage caused when the denotation is disregarded and improper use is performed.

|                                                                                                  |                                                                          |
|--------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------|
|  <b>DANGER</b>  | Denotes a potential hazard that will result in serious injury or death.  |
|  <b>WARNING</b> | Denotes a potential hazard that could result in serious injury or death. |
|  <b>CAUTION</b> | Denotes a hazard that could result in minor injury or property damage.   |

### **WARNING**

#### ► Preventing accidents

- Do not store within the reach of children or infants. Do not let them use it.
  - Doing so may cause an accident or injury due to accidental ingestion of removable parts.

#### ► USB power adaptor/USB cable

- Do not immerse them in water or wash with water.
- Do not place them over or near water filled sink or bathtub.
- Never use them if they are damaged or if the power plug fits loosely in a household outlet. Do not damage them.
- Do not deface, modify, forcefully bend or pull them. Also, do not place anything heavy on, twist or pinch them.
  - Doing so may cause electric shock or fire due to a short circuit.

- Do not use anything other than the supplied USB power adaptor or USB cable for any purpose. Also, do not use any other product with the USB power adaptor or USB cable. (See page 7.)
  - Doing so may cause burn or fire due to a short circuit.
- Do not connect or disconnect them to a household outlet with a wet hand.
  - Doing so may cause electric shock or injury.
- Regularly clean them to prevent dust from accumulating.
  - Failure to do so may cause fire due to insulation failure caused by humidity. Remove the USB power adaptor and wipe it with a dry cloth.
- Be sure to wipe off any water droplets from the socket before charging.
  - Failure to do so may cause electric shock, injury or fire.

#### ► In case of an abnormality or malfunction

- Immediately stop using and remove the USB power adaptor.
  - Failure to do so may cause fire, electric shock, or injury.
- <Abnormality or malfunction cases>
  - The main unit, USB power adaptor or USB cable is deformed or abnormally hot.
  - The main unit, USB power adaptor or USB cable smells of burning.
  - There is an abnormal sound during use or charging of the main unit.
  - Immediately request inspection or repair at an authorized service center.

#### ► Take care of the following precautions as well

- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.
  - Failure to do so may cause an accident or injury.

## **WARNING**

### ► **This product**

- This product has a built-in rechargeable battery. Do not throw into fire, apply heat, or charge, use, leave, or store in a high temperature environment.
  - Doing so may cause overheating, ignition, or explosion.
- Do not modify or repair.
  - Doing so may cause fire, electric shock, or injury.
- Contact an authorized service center for repair (replacement of battery etc).
- Never disassemble except when disposing of the product.
  - Doing so may cause fire, electric shock, or injury.

### ► **Power supply**

- Do not use in a way that would exceed the rating of the household outlet or the wiring.
  - Exceeding the rating by connecting too many plugs to one household outlet may cause fire due to overheating.
- Fully insert the power plug of USB power adaptor.
  - Failure to do so may cause fire or electric shock.
- Disconnect the USB plug when cleaning the main unit with water.
  - Failure to do so may cause electric shock or injury.

### ► **Cleaning**

- Always unplug the adaptor from a household outlet when cleaning the USB power adaptor and USB cable.
  - Failure to do so may cause electric shock or injury.

## **CAUTION**

### ► **Preventing damage to teeth and gums**

- People who are unable to operate should not use this appliance.
- People with weak oral senses should not use this appliance.
  - Doing so may cause injury or damage to teeth and gums.

- People who may have periodontal disease, have teeth being treated, or are concerned about symptoms within their mouth should consult a dentist before use.
  - Failure to do so may cause injury or damage to teeth and gums.

### ► **This product**

- Never use for a purpose other than oral cleaning.
  - Doing so may cause an accident or injury when the stream of water is pointed toward the eyes, nose, ears, tongue, or throat.
- Do not fill the tank with water warmer than 40 °C.
  - Doing so may cause burns.
- Do not share the nozzle with your family or other people.
  - Doing so may cause infection or inflammation.

### ► **Note the following precautions**

- Do not drop or subject to shock.
  - Doing so may cause injury.
- Do not allow metal objects or trash to adhere to power plug of USB power adaptor or USB plug.
  - Doing so may cause electric shock or fire due to a short circuit.
- Do not charge when the appliance socket is wet.
  - Doing so may cause electric shock or fire.
- Do not wrap the cord around the USB power adaptor when storing.
  - Doing so may cause the wire in the cord to break with the load, and may cause fire due to a short circuit.
- Disconnect the USB power adaptor or the USB plug (Type-C) by holding onto it instead of the cord.
  - Failure to do so may cause burns or electric shock due to breakage to the wire of the USB cable, or result in fire due to a short circuit.
- Disconnect the USB power adaptor from the household outlet when not charging.
  - Failure to do so may cause electric shock or fire due to electrical leakage resulting from insulation deterioration.
- After use, clean the appliance and the nozzle and dry them.
  - Failure to do so may cause damage to your health due to propagation of mould and bacteria.

## Handling of the removed battery when disposing

### **DANGER**

- The rechargeable battery is exclusively for use with this appliance.  
Do not use the battery with other products.
- Do not charge the battery after it has been removed from the product.
- Do not throw into fire or apply heat.
- Do not hit, disassemble, modify, or puncture with a nail.
- Do not let the positive and negative terminals of the battery get in contact with each other through metallic objects.
- Do not carry or store the battery together with metallic jewellery such as necklaces and hairpins.
- Do not use or leave the battery where it will be exposed to high temperatures, such as under direct sunlight or near other sources of heat.
- Never peel off the tube.  
- Doing so may cause overheating, ignition, or explosion.

### **WARNING**

- After removing the rechargeable battery, do not keep it within reach of children and infants.  
- The battery will harm the body if accidentally swallowed.  
If this happens, consult a physician immediately.
- If the battery fluid leaks out and comes in contact with the eyes, do not rub the eyes, and thoroughly rinse with clean water, such as tap water.  
- Failure to do so may cause eye injury.  
Consult a physician immediately.

### **CAUTION**


- If the battery fluid leaks out and comes in contact with the skin or clothes, rinse off with clean water, such as tap water.  
- Failure to do so may cause inflammation.

## Important information

- Please use the appliance together with toothbrush.  
It is hard to remove tartar or dental calculus attached to the tooth surface just by a water jet, so a toothbrush is required to brush the teeth. We recommend rinsing the inside of your mouth with the water jet after brushing your teeth.
- When you first begin using the oral irrigator, you may experience slight bleeding of the gums even though your gums may be healthy. This is because your gums are receiving stimulation by the oral irrigator for the first time. The bleeding should stop after one or two weeks of use. Bleeding that persists for more than two weeks may be indication of an abnormal condition of your gums, such as periodontal disease. In this case, stop using the oral irrigator and consult your dentist.
- Panasonic's specified adaptor with splash-resistant design and a USB cable with an overheating protection feature should be used to ensure safety. The adaptor is not completely waterproof, so prevent it from getting wet.
- Do not put salt, mouthwash or other chemicals into the water tank.  
Doing so may cause malfunction.
- Do not dry the main unit with a dryer or a fan heater.  
Doing so may cause malfunction or deformation of the parts.
- Oral irrigators shall be used only with cold or warm water, or specific solutions as defined in the manufacturer's instructions.
- Children shall not play with the appliance.
- Store the appliance in a place with low humidity after use. Leaving the appliance in a bathroom may cause malfunction.
- Charge the appliance correctly according to these Operating Instructions. (See pages 8 and 9 "Charging".)
- This appliance contains a battery that is only replaceable by skilled persons. Please contact an authorized service center for the details of repair.
- The following symbol indicates that a specific detachable power supply unit is required for connecting the electrical appliance to the supply mains. The type reference of power supply unit is marked near the symbol.

## Parts identification
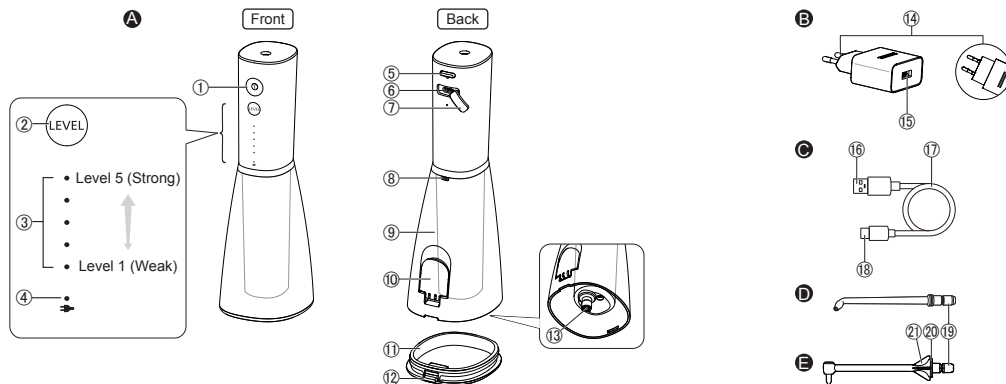
### A Main unit

- ① Power Switch
- ② Level switch
- ③ Water pressure level indicator (blue)
- ④ Charge status indicator (red)
- ⑤ Nozzle release button
- ⑥ Socket
- ⑦ Socket cover
- ⑧ Air vent
- ⑨ Water tank (inside)
- ⑩ Water tank cap
- ⑪ Bottom cover
- ⑫ Bottom cover release button
- ⑬ Filter

### B USB power adaptor (ES-RU1-02)

- ⑭ Power plug
- ⑮ USB port (Type-A)

### C USB cable

- ⑯ USB plug (Type-A)
- ⑰ Cord
- ⑱ USB plug (Type-C)

### D Ultrasonic nozzle (×1)

### E Orthodontic nozzle (×1)

- ⑲ Nozzle insertion port
- ⑳ Identification ring
- ㉑ Nozzle Handle
- The nozzle is a consumable.  
(See page 13.)

### Note

There may be moisture inside the main unit (water tank) left from the distilled water used for product testing.


## Charging

- Please wipe off any water droplets from the socket or socket cover.

### 1 Open the socket cover, and insert the USB plug (Type-C) of the USB cable into the socket on the main unit.

#### ► If the socket cover comes off


- 1 While holding the main unit with your hand, press and hold the socket cover with your thumb so that the soft rubber joint (mushroom-shaped) of the socket cover is aligned with the fixing hole.
- 2 As shown in the figure, push the soft rubber joint (mushroom-shaped) into the fixing hole with a SIM ejector pin. (Do not use a sharp-pointed tool as this may cause damage to the soft rubber joint or a hand injury.)
- 3 Make sure that the soft rubber joint (mushroom-shaped) is fully pushed into the fixing hole.
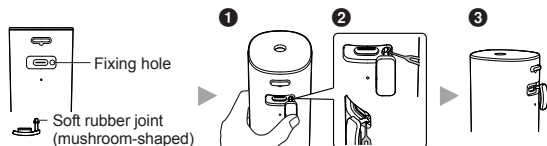
### 2 Insert the USB plug (Type-A) of the USB cable into the USB port (Type-A) of the USB power adaptor.

- Check the direction of the USB plug (Type-A) terminals, hold the USB plug (Type-A), and insert or pull it out straight. (Deflection or insertion in a wrong direction may cause deformation or malfunction.)

### 3 Insert the USB power adaptor into the household outlet.
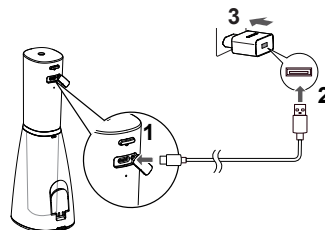
### 4 Charging is completed after approx. 1 hour. For safety, disconnect the USB power adaptor from the household outlet after charging is completed. (The charging time may vary with the remaining battery capacity.)

- Remove the USB plug (Type-C) from the socket, and close the socket cover.
- Charging the battery for more than 1 hour will not affect battery performance.
- This appliance uses a Lithium-ion battery. Therefore, the service life of the battery will not be affected even if it is charged after each use.

### ► Operating time after a full charge

**Approx. 10 minutes\***

(When the water pressure level is set to 5)

\* It varies with the ambient temperature and usage conditions.

### ► About the charging reminder function

- If the battery is running low, the charge status indicator will start flashing (once a second) during use and light off in 5 seconds after turning off the power switch.

### ► Display condition of the charge status indicator

|                      |                                              |
|----------------------|----------------------------------------------|
| Charging in progress | Lights on                                    |
| Charging complete    | Lights off                                   |
| Charging abnormal    | Flashing<br>(Flashing in red twice a second) |

### ► How to check if charging is completed

- After charging is completed, if you disconnect and then connect the power supply, the charge status indicator will light on and then light off in about 5 seconds.


### ► Notes

- The charging time may differ when the appliance is used for the first time, or if it has not been used for more than 6 months. In addition, the charge status indicator may not light on in the first several minutes of charging, but it will light on as the charging proceeds.
- The recommended ambient temperature for charging is 10 °C - 35 °C. It may take longer to charge, or the charge status indicator will flash rapidly and the battery may fail to charge properly in extremely high or low temperatures.
- The battery will deteriorate if not used for over 6 months, so make sure to perform a full charge at least once every 6 months.

## How to use


### 1 Fill the water tank with water.

- 1 Open the water tank cap.
- 2 Add water while keeping the main unit level.
  - Do not fill beyond the full-water position.
- 3 Close the water tank cap securely.  
(When the full-water position is reached, or when the water tank cap is closed after the full-water position is reached, a small amount of water may flow out of the air vent.)

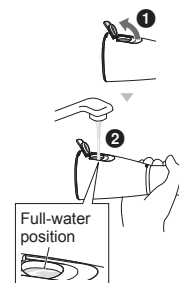
#### ► If the water tank cap comes off

You can attach the water tank cap by following the procedure below when it comes off.

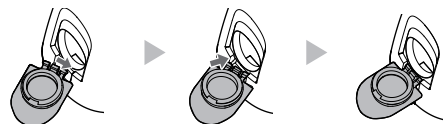
### 2 Install the nozzle.

- Choose a nozzle for your purpose.  
(See page 11.)
- Press in firmly until a "click" is heard.
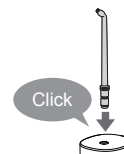
### 3 Press the level switch, and select the water pressure level.

- Each time you press the level switch, the water pressure level is switched as follows:  
1 (weak) → 2 → 3 → 4 → 5 (strong) → 1...

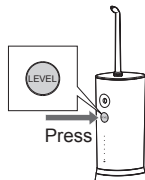► **About the water pressure level memory function**  
The last selected water pressure level will be remembered, from which the appliance will start next time.

- When the appliance has not been used for a long time or when the battery is used up, the water pressure level will return to 1.

### 4 Place the nozzle in your mouth, close your mouth gently, and then turn on the power switch to rinse your teeth.

- You can switch the water pressure level by pressing the level switch, even while using the appliance.
- If the main unit is tilted too much, the water does not come out from the nozzle.
- During use, be careful not to cover the air vent with your hand. Otherwise, the water does not come out from the nozzle.

Keep your face over the basin.

\*It can also be used in a bath.

► **Operation time with full water approx. 60 seconds.**

(When the water pressure level is set to 5)

### 5 To prevent the water from splashing about, turn off the power switch while keeping the nozzle in your mouth.

- Please clean the appliance after each use. (See page 12.)

#### ► Notes

- The appliance does not work during charging.
- Always use new tap water or lukewarm water. (Do not put salt, mouthwash or other chemicals into the water tank.)
- If you feel the water is cold, use lukewarm water with a temperature of 40 °C or less.
- For the first use, turn on the power switch to jet water so that you can check the water pressure level. (Please start using the appliance from water pressure level 1.)
- The appliance may not operate in an ambient temperature approximately 5 °C or lower.
- Do not push the tip of the nozzle too hard against the teeth or gums. (Except the orthodontic nozzle)
- When you feel that the irritation is strong, turn off the power switch and use the appliance from water pressure level 1.
- Do not turn on the power switch when the water tank is empty unless you are drying the moisture inside the main unit. (Doing so may cause malfunction.)

## Points of care

- Adjust the angle and movement of the nozzle according to the purpose and the area you want to clean.
- For areas where the gums are sensitive, please set the water pressure level to 1.

### Ultrasonic nozzle

To clean with powerful water jet.

### How to use

#### Interdental cleaning

- Direct the water jet onto the space between the teeth.
- Direct from the front and rear sides of the teeth.

#### Periodontal pocket cleaning

- Aim the water jet at the space between teeth and gums.
- Slowly shift the water jet along the gums.
- Clean the front and back of the teeth and behind the molars as well.

#### Gums care

- Stimulate the gums.

### Orthodontic nozzle

To clean with 2-in-1 water jet and brush nozzle.

### How to use


#### Lightly place the brush tip on the crooked teeth.

- Place the brush tip on the area that is difficult for a toothbrush to reach and tends to remain unpolished.
- Wiggle the brush mainly in interdental spaces.

#### To adjust the orientation of the nozzle


You can adjust the direction of water jet by rotating the nozzle handle.
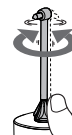
#### If the identification ring comes off

Align the nozzle with the identification ring (in any direction) to install it.

- Even if the identification ring is not attached, the performance of the nozzle is not affected.

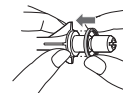
### For orthodontic bracket area -

#### Ultrasonic nozzle

- Direct the water jet around orthodontic brackets and bridges.
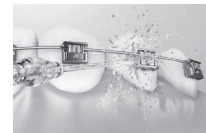
#### Orthodontic nozzle

- Insert the brush to the gap between the teeth and the orthodontic bracket.
- Wiggle the brush tip.



## Cleaning

- Please clean the appliance after each use.

### 1 Pour out the remaining water.


- 1 Open the water tank cap, and pour out the remaining water.
  - Be sure to pour out the remaining water after each use.
- 2 Turn on the power switch to drain the remaining water in the main unit.
- 3 Turn off the power switch when no more water comes out.
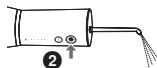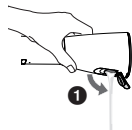
### 2 Remove the nozzle and bottom cover. Wash the nozzle, main unit and inside of the tank with water.

#### ► How to remove the nozzle


Pull out the nozzle while pressing the nozzle release button.
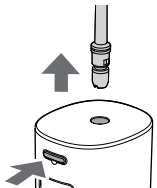
#### ► How to remove the bottom cover

Bottom cover release button

Lift up the bottom cover while pressing the button.

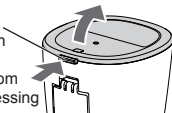
### Preparations

- Make sure the socket cover is closed before washing with water.

### Washing

- Do not hit the product on the washstand when draining.

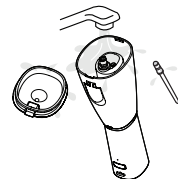
### 3 Store the appliance after drying.

- Wipe away the water with a towel and dry them in a well-ventilated shade.
- Attach the bottom cover after it is air-dried.

#### ► How to attach the bottom cover

- 1 Insert the bayonet on the front end of the bottom cover into the main unit.
- 2 Close the bottom cover while pressing the bottom cover release button.
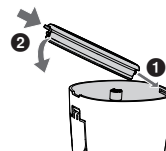
#### ► If the rubber ring of the bottom cover comes off

Install the rubber ring into the groove of the bottom cover.


### ► Notes

- Do not wash with lukewarm water exceeding 50 °C.
- Use neutral kitchen detergent and rinse well with water to ensure no detergent remains after cleaning.
- Never use thinners, benzine, or alcohol.  
(Doing so may cause malfunction, or cracking or discolouration of the parts.)

### ► When the water jet is poor

Remove the dirt from the filter.

Filter
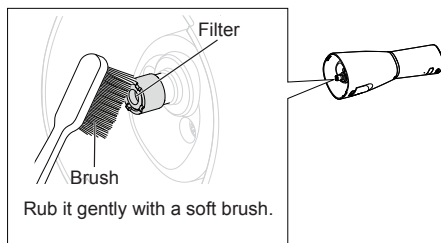
### ► When dirt is a concern

USB power adaptor and USB cable

Wipe it off with a dry cloth. (Do not wash with water)

## Changing the nozzle (Sold separately)

- The nozzle is a consumable.

| Part number | Replacement Nozzle (×2) |  |
|-------------|-------------------------|--|
| WEW0986     | Ultrasonic nozzle       |  |
| WEW0984     | Orthodontic nozzle      |  |

### ■ Ultrasonic nozzle

- Guidelines for replacement: approximately 6 months.  
(For hygiene reasons)
- Once the tip of the nozzle gets deformed, replace the nozzle with a new one even within 6 months.

### ■ Orthodontic nozzle

- Guidelines for replacement: approximately 3 months.
- Once the brush bristles are bent out, replace the nozzle with a new one even within 3 months.
- The replacement nozzle comes with an identification ring so that a same product can be shared among several people.

|                    | New nozzle (tip) | Deformed (worn) nozzle (tip) |
|--------------------|------------------|------------------------------|
| Ultrasonic nozzle  |                  |                              |
| Orthodontic nozzle |                  |                              |

## Troubleshooting

|                | Problem                                                                                     | Possible cause                                                                                                               | Action                                                                                                                                                                                                                                                 |
|----------------|---------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| White charging | The charge status indicator does not light on.                                              | ▶ The appliance was just purchased or it has not been used for a while.                                                      | ▶ Keep charging, and the indicator will light on after a while. If it does not light on after a while, please request a repair at your dealer or Service Center. (See pages 8-9.)                                                                      |
|                |                                                                                             | ▶ When the socket is wet or foreign objects adhered to the socket, the USB cable will activate the self-protection function. | ▶ Disconnect the power supply, check the socket, and remove any water or foreign objects from the socket.                                                                                                                                              |
|                | The charge status indicator flashes rapidly (twice a second).                               | ▶ The USB cable is not firmly connected to the USB power adaptor.                                                            | ▶ Check to see that the USB cable is in full contact with the USB power adaptor, and that the charge status indicator is on. (See page 9.)                                                                                                             |
|                |                                                                                             | ▶ The ambient temperature is too low or too high during charging.                                                            | ▶ Please charge within the recommended temperature range of 10 °C to 35 °C.<br>▶ If the indicator flashes even when you charge it within the recommended temperature range, contact your dealer or Service Center for repair.                          |
|                | The appliance can be used for just several minutes even after a full charge.                | ▶ The USB cable is not firmly connected to the USB power adaptor.                                                            | ▶ Make sure that the USB plug (Type-C) of the USB cable is fully inserted into the socket of the main unit or the USB plug (Type-A) of the USB cable is fully inserted into the power supply and the charging status indicator is on. (See pages 8-9.) |
|                |                                                                                             | ▶ The battery has reached the end of its life (approx. 3 years).                                                             | ▶ Replace it with a new one.<br>▶ Please contact your dealer or Service Center.                                                                                                                                                                        |
|                | The appliance cannot be charged normally, but the charging status indicator does not flash. |                                                                                                                              |                                                                                                                                                                                                                                                        |

11/11/2016

## Battery life

Approximately 3 years. The battery may have reached its operating life if the usage frequency decreases drastically, even after a full charge. (The operating life of the battery may differ depending on usage or storage conditions.)

## Removing the built-in rechargeable battery


**Remove the built-in rechargeable battery before disposing of the appliance.**

The battery is to be disposed of safely.

Please make sure that the battery is disposed of at an officially designated location if there is one.

This figure must only be used when disposing of the appliance, and must not be used to repair it. If you dismantle the appliance yourself, it will no longer be waterproof, which may cause it to malfunction.

- Remove the appliance from the USB cable and USB power adaptor, press the power switch to turn on the power and then keep the power on until the battery is completely discharged.
- Remove the bottom cover.
- Perform steps ① to ⑥ to disassemble the main unit and remove the battery with a screwdriver.
- Insulate the ⊕ and ⊖ terminals by applying tape to them. (To prevent short circuit of terminals).
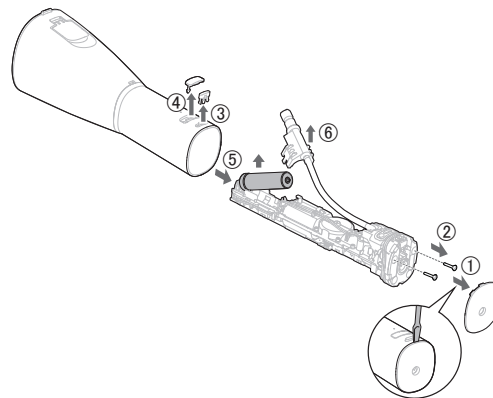
### For environmental protection and recycling of materials

This appliance contains a Lithium-ion battery.

Please make sure that the battery is disposed of at an officially assigned location, if there is one in your country.

## Specifications

|                |                                                                                                      |
|----------------|------------------------------------------------------------------------------------------------------|
| Power source   | See the name plate on the product<br>(Automatic voltage conversion)                                  |
| Charging time  | Approx. 1 hour                                                                                       |
| Operating time | LEVEL 5: Approx. 10 minutes*<br>(*at 20 °C when fully charged)                                       |
| Tank volume    | Approx. 200 mL<br>(For approx. 60 seconds of operation when the<br>water pressure level is set to 5) |
| Dimensions     | Main unit: 212 (H) x 68 (W) x 81 (D) mm<br>(Not including nozzle)                                    |
| Mass           | Main unit: Approx. 275 g (Including nozzle)                                                          |

This product is intended for household use only.

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## This image shows a blank sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

## This image shows a blank sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

Panasonic Corporation  
<https://www.panasonic.com>  
© Panasonic Corporation 2024

© EN, VT

Printed in China  
EW9700DJ861442 S0824-0