Panasonic

Operating Instructions (Household) IPL Hair Removal System Model No. ES-WH80

English..... 3

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Operating Instructions (Household) IPL Hair Removal System Model No. ES-WH80

Contents

Safety precautions.....	4	For healthy skin	14
Intended use	9	Lamp indication	14
Parts identification	9	Cleaning	16
Treatment areas	10	Lamp cartridge replacement	16
Guide for treatment frequency, time..	10	Q&A.....	17
Charging.....	11	Troubleshooting	18
Before use	12	Battery life	19
When using the appliance for the first time	12	Removing the built-in rechargeable batteries.....	19
How to use	12	Specifications	20

Thank you for purchasing this Panasonic product.

Before operating this unit, please read these instructions completely and save them for future use.

Safety precautions

To reduce the risk of injury, loss of life, electric shock, fire, and damage to property, always observe the following safety precautions.

Explanation of symbols

The following symbols are used to classify and describe the level of hazard, injury, and property damage caused when the denotation is disregarded and improper use is performed.

 DANGER	Denotes a potential hazard that will result in serious injury or death.
 WARNING	Denotes a potential hazard that could result in serious injury or death.
 CAUTION	Denotes a hazard that could result in minor injury.

The following symbols are used to classify and describe the type of instructions to be observed.

	This symbol is used to alert users to a specific operating procedure that must not be performed.
	This symbol is used to alert users to a specific operating procedure that must be followed in order to operate the unit safely.

WARNING

► To prevent skin or health problems

 **Persons using an implanted medical electronic device, such as a pacemaker, should not allow this unit to touch the skin where the medical electronic device is implanted.**

- Doing so may cause an accident or illness.

 **Do not use during pregnancy or menstruation.**
- Doing so may cause skin problems due to unstable hormone balance.

 **Do not use on tattoos.**
- Doing so may cause keloids by burns or redness or other rashes on skin.

 **Do not shine directly into the eyes.**
- Doing so may cause injury to the eyes.

 **Only use on underarms, arms, legs or bikini-line.**
- Failure to do so may cause burns or redness or other rashes on skin.

► Power supply

 **Do not connect or disconnect the adaptor to a household outlet with a wet hand.**
- Doing so may cause electric shock or injury.

WARNING

Do not damage or modify, or forcefully bend, pull, or twist the cord. Also, do not place anything heavy on or pinch the cord.

Do not immerse the AC adaptor in water or wash it with water.

 **Do not place the AC adaptor over or near water filled sink or bathtub.**

Never use the appliance if the AC adaptor is damaged or if the power plug fits loosely in a household outlet.

- Doing so may cause electric shock or fire due to a short circuit.

 **Do not use anything other than the supplied AC adaptor. Also, do not charge any other product with the supplied AC adaptor.**

- Doing so may cause burn or fire due to a short circuit.

 **Do not use in a way exceeding the rating of the household outlet or the wiring.**

- Exceeding the rating by connecting too many plugs to one household outlet may cause fire due to overheating.

 **Always ensure the appliance is operated on an electric power source matched to the rated voltage indicated on the AC adaptor.**

Fully insert the adaptor or appliance plug.

- Failure to do so may cause fire or electric shock.

 **Always unplug the adaptor from a household outlet when cleaning.**

- Failure to do so may cause electric shock or injury.

Regularly clean the power plug and the appliance plug to prevent dust from accumulating.

 - Failure to do so may cause fire due to insulation failure caused by humidity.

Disconnect the adaptor and wipe with a dry cloth.

► This product

 **This product has a built-in rechargeable battery. Do not throw into flame or heat. Do not charge, use, or leave in high temperatures.**

- Doing so may cause overheating, ignition, or explosion.

 **Do not store within the reach of children. Do not let them use it.**

- Doing so may cause an accident or injury.

Do not modify or repair.

 - Doing so may cause fire, electric shock, or injury.
Contact an authorized service center for repair (battery change etc).

 **Never disassemble except when disposing of the product.**

- Doing so may cause fire, electric shock, or injury.

 **Do not store the appliance in a bathroom or location with high humidity. Also, do not place in a location where water may spill (such as top of a vanity).**

- Doing so may cause electric shock or fire.

Do not immerse appliance in water or wash it with water.

 - Doing so may cause electric shock or fire due to a short circuit.

Keep the appliance dry.

! WARNING

► In case of an abnormality or malfunction

Immediately stop using and remove the adaptor if there is an abnormality or malfunction.

- Using it in such conditions may cause fire, electric shock, or injury.

<Abnormality or malfunction cases>

- ! • The main unit, adaptor or cord is deformed or abnormally hot.
- The main unit, adaptor or cord smells of burning.
- There is abnormal sound during use or charging of the main unit, adaptor or cord.
- Immediately request inspection or repair at an authorized service center.

► Take care of the following precautions as well

Do not remove the lamp cartridge except during replacement.

Do not touch the  area when replacing the lamp cartridge.

- ⊘ - Doing so may cause electric shock or fire due to a short circuit.

This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

Children should be supervised to ensure that they do not play with the appliance.

- Failure to do so may cause an accident or injury.

The supply cord cannot be replaced. If the cord is damaged, the AC adaptor should be scrapped.

- Failure to do so may cause an accident or injury.

CAUTION

► To protect the skin

The following persons should not use this appliance:

- Persons with dermatitis
- Persons suffering from atopic dermatitis or rashes and other skin diseases
- Persons with allergic diathesis, or sensitive skin which is easily irritated by cosmetics, clothing, metals, etc.
- Persons with skin that is sensitive to sunlight

Do not use in the following body area:

- On suntanned skin
- Dark brown or black skin areas such as moles (Refer to the sheet for "Suitable Hair and Skin Tones".)
- Scabs, injuries, rashes, bruises, warts, pimples, skin cancer, ringworm (dermatophytosis), folliculitis, irritated skin, insect bites, etc.
- Face and genitals
- Angioma
- Areas with cosmetic surgery

Do not use in the following cases:

- When you are not well or when you feel an abnormality in your skin or body

Do not use on the same area repeatedly, or use IPL hair removal more than 3 times per week.

- Doing so may cause injury to the skin, inflammation, bleeding, increased pain or other problems or aggravated symptoms.

Do not use a cracked lamp cartridge.

- Doing so may cause injury to the skin.

Do not remove the frame except during cleaning or lamp cartridge replacement.

- Doing so may cause injury to the skin or fingers.

Only for use by persons over the age of 18.

- Failure to do so may cause injury or inflammation of the skin, bleeding or increased pain.

Before use, check that there is no distortion, cracks or damage to the frame and flash window.

- Failure to do so may cause injury to the skin.

Persons on medication should consult a dermatologist.

-  - Skin may become sensitive to light stimulation depending on the type of medication being taken, and may cause skin problems.

If inflammation (such as skin turning red) occurs, stop using the appliance, and consult a dermatologist if the inflammation continues even after 2 days.

- Failure to do so may cause symptoms to worsen.

► Power supply

Do not allow pins or trash to adhere to the power plug or the appliance plug.

-  - Doing so may cause electric shock or fire due to a short circuit.

Do not wrap the cord around the adaptor when storing.

-  - Doing so may cause the wire in the cord to break with the load, and may cause fire due to a short circuit.

Disconnect the adaptor or the appliance plug by holding onto the adaptor or the appliance plug instead of the cord.

-  - Disconnecting by holding onto the cord may cause electric shock or injury.

CAUTION

-  **Disconnect the adaptor from the household outlet when not charging.**
- Failure to do so may cause electric shock or fire due to electrical leakage resulting from insulation deterioration.

► This product

-  **Do not use on skin that is wet with water, lotion, shaving cream or other products.**
- Doing so may cause an electric shock or malfunction.

-  **Do not drop or subject to shock.**
- Doing so may cause injury.

-  **Do not remove the frame immediately after using IPL hair removal.**
- Doing so may cause burns as the flash window becomes hot.

► Note the following precautions

-  **Do not use on clothes or head hair.**
- Doing so may cause burns or discoloration.

► Disposing of the rechargeable battery

DANGER

The rechargeable battery is exclusively for use with this appliance. Do not use the battery with other products. Do not charge the battery after it has been removed from the product.

- Do not throw into fire or apply heat.
- Do not hit, disassemble, modify, or puncture with a nail.
-  • Do not let the positive and negative terminals of the battery get in contact with each other through metallic objects.
- Do not carry or store the battery together with metallic jewellery such as necklaces and hairpins.
- Do not charge, use, or leave the battery where it will be exposed to high temperatures, such as under direct sunlight or near other sources of heat.
- **Never peel off the tube.**
- Doing so may cause overheating, ignition, or explosion.

WARNING

After removing the rechargeable battery, do not keep it within reach of children and infants.

-  - The battery will harm the body if accidentally swallowed. If this happens, consult a physician immediately.

! WARNING


If the battery fluid leaks out, take the following procedures. Do not touch the battery with your bare hands.

- The battery fluid may cause blindness if it comes in contact with your eyes.
Do not rub your eyes. Wash immediately with clean water and consult a physician.
- The battery fluid may cause inflammation or injury if it comes in contact with the skin or clothes.
Wash it off thoroughly with clean water and consult a physician.

Intended use

- Use IPL hair removal only after removing excess hair. Emitting light without removing hair may cause hair to become trapped in the guard or flash window, and the dirt cannot be removed. Before use, wipe away any water or sweat from the treatment area, and any hair sticking to the skin first.
- Do not use thinner, benzene, alcohol, etc. to clean the appliance. Doing so may cause failure, cracking or discoloration of the main body. Wipe with a dry cloth.
- Do not leave the appliance in an area with high temperature. Doing so may cause malfunction.

Parts identification

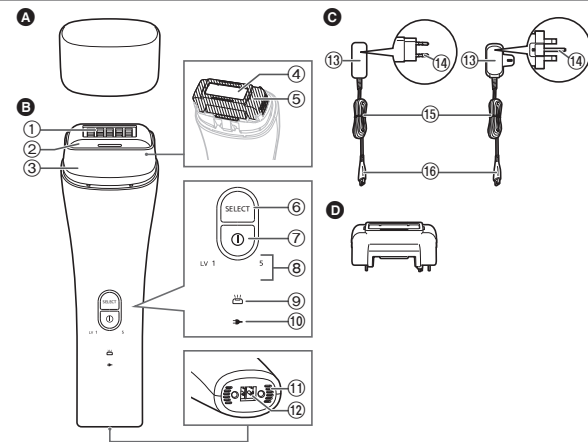A Cap

B Main body

- Guard (part with grid pattern)
- Flash switch
- Frame
- Flash window
- Lamp cartridge
- Level selector switch (flashes)
- Power switch
- Level indication lamp (blue) (flashes) (See pages 12 and 14.)
- Cartridge confirmation lamp (red) (See page 14.)

- Charging indication/notification lamp (red) (See page 14.)

- Air inlet
- Socket

C AC adaptor (RE7-87) (The shape of the AC adaptor differs depending on the area.)

- Adaptor
- Power plug
- Cord
- Appliance plug

Accessory


- Lamp cartridge (x1)

Treatment areas

The appliance can be used on the  areas.

Underarms

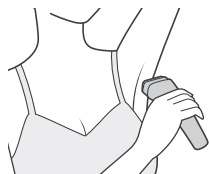
Arms

Legs

Bikini-line*

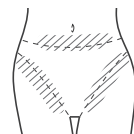* For bikini-line treatment, wear a swim wear, under wear etc. first to check the treatment area.

Guide for treatment frequency, time

Treatment frequency

Legs/Arms

Use once every 2 weeks for the first 2 months.

Use once every 4 weeks from the 3rd month onward.

Underarms/Bikini-line

Use once every week for the first 2 months.

Use once every 2 weeks from the 3rd month onward.

Guide for treatment time

Legs	Approximately 10 minutes per session for both legs (approximately 190 flashes)
Arms	Approximately 5 minutes per session for both arms (approximately 90 flashes)
Underarms	Approximately 2 minutes per session for both underarms (approximately 30 flashes)
Bikini-line	Approximately 2 minutes per session (approximately 30 flashes)

- Do not use on the same area repeatedly, or use IPL hair removal more than 3 times per week.
Try using the appliance up to a maximum of 2 times per week even if removal of unwanted hair is not immediately obvious. (Treatment effects may differ depending on the area or individual body types.)
- If the appliance has not been used for a while, resume use at the same frequency as the first 2 months.
- Treatment time may differ depending on the area or amount of hair, treatment area and how the appliance is used.
- Please remove any unwanted hair first from the treatment area. Please use a separate razor or electric shaver.


- The appliance differs from permanent hair removal treatment. Unwanted hair will gradually grow back if treatment is stopped.

Charging

- 1** Insert the appliance plug firmly into the socket at the bottom of the main body.
- 2** Firmly insert the adaptor into a household outlet.
- 3** Check that the charging indication/notification lamp glows.

While charging	After charging is completed
	
Red light glows.	Red light goes off.

- The charging time may differ when the appliance is being used for the first time, or if it has not been used for half a year or more. The lamp on the main body may not glow for several minutes after charging starts, but will glow if charging continues.
- Charging is completed in about 3 hours.
- The charging time may differ with the charging capacity.
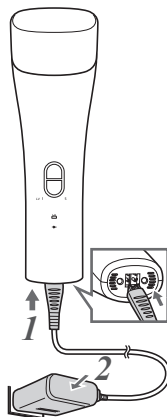
<To check if charging is completed>
Press the power switch with the AC adaptor connected. The charging indication/notification lamp glows, and goes off after approximately 5 seconds.

- 4** Remove the adaptor from the household outlet, and remove the appliance plug from the main body.

(for safety and to reduce energy usage)

- There will be no problems with the main body even if it is charged for 3 hours or more.
- The appliance can be used for approximately 600 flashes with a full single charge.
- The number of flashes may differ depending on the ambient temperature and how the appliance is used.

Notes

- Recommended ambient temperature for charging is 15 – 35 °C. Charging may not be possible outside of the recommended temperature.
- Charging the appliance every time it is used will not affect the life of the rechargeable battery.
- If there is crackling noise on radios or other devices while charging, charge the appliance using a different household outlet.

Before use

- Do not use the appliance in dark places. Doing so may cause temporary irritation in the eyes.
- Please remove any unwanted hair first from the treatment area.
- Before use, wipe away any water or sweat from the treatment area, and any hair sticking to the skin. Otherwise, the dirt on the guard or flash window cannot be removed.

When using the appliance for the first time

Check that the output level suits your skin.

Set the output level to 1 (LV1) and flash the treatment area once.

After 2 days

If there are no skin abnormalities

Set the output level to 1 (LV1) and start treatment on the area that you want.

- Guide for treatment frequency, time (See page 10.)

After 2 weeks*


If there are no skin abnormalities

*For armpits and bikini line: after 1 week

To increase the output level for treatment, flash once on the treatment area, and check your skin condition for around 2 days before treatment. Increase the output level one step at a time.

- If there is intense pain or redness after 2 or more days, stop using the appliance and consult a dermatologist.

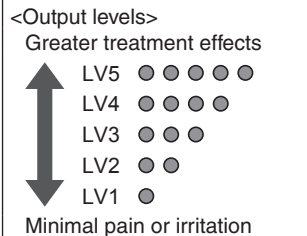How to use

- The power will not turn on even when pressing the switch with the AC adaptor connected to the main body.
- Recommended ambient temperature for operation is 5 – 35 °C. The appliance might not be usable outside of the recommended temperature.

1 Remove the cap.

2 Press and hold (for approximately 3 seconds) to turn the power on. ①

- The power will turn off automatically if the appliance is not used for approximately 1 minute.


3 Press  to select the output level.

- Level 1 is selected when the power is turned on.
- The level cycles through 1 → 2 → 3 → 4 → 5 → 1...each time the level selector switch is pressed.
- The level indication lamp blinks, and then turns on. (A beep sounds, indicating that preparation for flash is complete.)

4 When the level indication lamp glows, apply the guard closely against the skin and keep it pressed until flashes.

When the appliance beeps, press the head perpendicularly against the skin.

- Use the appliance after checking that there are no other people nearby.
- Points for IPL hair removal (See this page.)
- When the appliance beeps, preparation for flash is complete.

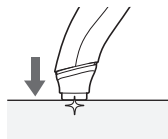
5 Press to turn the power off.

Points for proper IPL hair removal

Point 1


Apply the guard perpendicularly against the skin and keep it pressed against the skin.

The flash window will not flash unless it is applied perpendicularly to the skin.

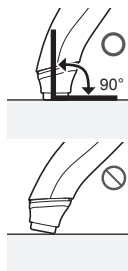
Point 2

The size of the guard (part with grid pattern) is a guide for the IPL hair removal range of each flash.

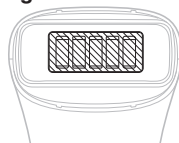
 :IPL hair removal range of each flash

- When the appliance beeps, preparation for flash is complete.

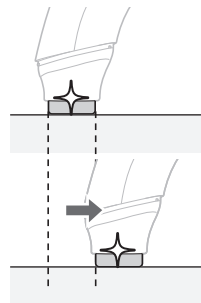Continuous flashes

After one flash, slide the appliance to the next flash area while it is still pressed against the skin. Up to 10 continuous flashes at an interval of approximately 3 seconds is possible.

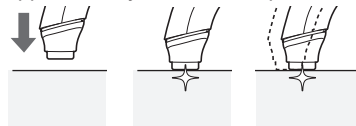- The flash window will no longer flash after 10 flashes. Remove the appliance from the skin and try again.
- Do not use the appliance continuously on the same area.

For healthy skin

The following care is recommended after IPL hair removal.

1. Cool the skin using a towel soaked in ice water.

2. Moisturize the skin with the lotion etc.

Lamp indication

► Normal indication

  	The level indication lamp (blue) blinks (once/second). The image on the left is when level 5 is selected. Preparation for flash in progress.
  	The level indication lamp (blue) glows. The image on the left is when level 5 is selected. Flashing is possible.

► Less than 1 000 flashes remaining until lamp cartridge replacement period

LV 1 5  	Cartridge confirmation lamp (red) glows. There are less than 1 000 flashes remaining until the lamp cartridge replacement period. Prepare a new lamp cartridge. (See pages 16 and 17.)
--	--

► Lamp cartridge replacement period


LV 1 5  	The cartridge confirmation lamp (red) blinks (once/second), and the power turns off 5 seconds later. IPL hair removal is not possible as the lamp cartridge replacement period is up. Replace the lamp cartridge. (See pages 16 and 17.)
--	--

► Charging notification

LV 1 5  	The charging indication/notification lamp (red) blinks (once/second). The battery level is low. Charging is recommended.
--	--

► Malfunction display

LV 1 5 	[Abnormal charging] The charging indication/notification lamp (red) blinks quickly (twice/second). Recommended ambient temperature for charging is 15 – 35 °C. Check the room temperature of the place where the appliance is going to be charged, and charge it again.
LV 1 5 	[Cartridge] The cartridge confirmation lamp (red) blinks quickly (twice/second), the appliance emits a series of beeps, and the power turns off after 5 seconds. The lamp cartridge is not attached correctly. Attach the lamp cartridge correctly after checking that the appliance has been turned off. (See page 16.)

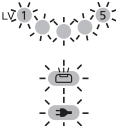	[Other malfunctions] The level indication lamp (blue), cartridge confirmation lamp (red) and charging indication/notification lamp (red) blink simultaneously (once/second), the appliance beeps twice, and the power turns off after 5 seconds. Recommended ambient temperature for operation is 5 – 35 °C. Check the room temperature of the place that the appliance is going to be used, and turn it on again. (The appliance may not adapt to the room temperature immediately after changing the location.) If the symptom occurs repeatedly, contact an authorized service center.
---	--


Cleaning

Always turn the main body power off and wait for 2 to 3 minutes before cleaning.

► Flash window/Frame

Remove any trimmings on the guard or flash window. Failure to do so will make it difficult to remove the dirt, and will weaken the light.

1. Remove the frame.

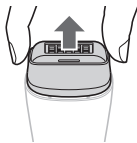2. Wipe any dirt clean carefully with a dry cloth.

Flash window

Frame/Guard

3. Attach the frame.
Insert until it clicks.

Lamp cartridge replacement

A guide for the lamp cartridge replacement is approximately 2 years when the appliance is used at the recommended frequency and time.

The lamp cartridge replacement period is approximately 20 000 flashes.


Replacement period notification

The cartridge confirmation lamp (red) glows when the power is turned on with less than approximately 1 000 flashes remaining. The cartridge confirmation lamp (red) blinks when turning the power on if the lamp cartridge replacement interval has come. (See page 14.)


Replacement and resetting of the lamp cartridge

1. Press  to turn the power off.
2. Remove the lamp cartridge.

- ① Remove the frame.

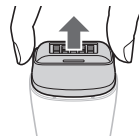- ② Remove it by pulling it straight in the direction of the arrow.

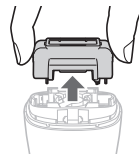
- Do not use a lamp cartridge that has been used. (Doing so may cause damage to the lamp cartridge.)

3. Attach a new lamp cartridge.

- ① Insert it straight in the direction of the arrow.

- Insert it firmly all the way in.

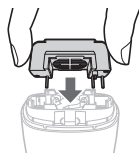- ② Attach the frame.

Insert until it clicks.

4. Press and hold (for approximately 3 seconds)  until the appliance beeps to reset it.

- Replacement is complete if the cartridge confirmation lamp (red) blinks twice and the alarm beeps.
- The power will not turn on unless the reset is completed.

Notes

- Always reset the appliance immediately after replacing the lamp cartridge. (The appliance cannot be reset if too much time has passed since replacing the lamp cartridge.)
- If the reset operation cannot be completed, recharge the appliance, reattach the lamp cartridge again and reset the appliance.

► Replacement parts

Replacement parts are available at your dealer or Service Center.

Replacement lamp cartridge	WES2W13
----------------------------	---------

Q&A

Question	Answer
What type of light is emitted from the flash window?	► Instantaneous flash from a xenon lamp.
Why can the appliance be used only by those over the age of 18?	► This is because of unstable hormone balance when growing up and the skin can become sensitive easily.
Can the appliance be used in areas other than underarms, arms, legs or bikini-line?	► Only use the appliance on underarms, arms, legs or bikini-line. (Do not use the appliance on the face or other areas.)
Is it alright to use the appliance everyday?	► Using the appliance on the same area everyday is not recommended. ► Please follow the guide for treatment frequency and time. (See page 10.)
Why is there a burning smell after a flash?	► This is because the flash applies heat to unwanted hair, and is not a sign of a malfunction. Remove any hair if it has grown too long. Also remove all trimmings remaining on the skin.
The light feels very bright while using the appliance but is this a problem?	► Apply the flash switch closely against the skin. (See page 13.) ► If the light is still a concern, wear sunglasses or other eye protection.

Question	Answer
The light feels hot while using the appliance but is this a problem?	▶ The sensation is different for each person but it may feel hotter if hair has not been removed properly. Please shave any hair short before flashing.
The treated skin turns red while using the appliance but is this a problem?	▶ The treated skin may temporarily turn red or feel painful depending on the body type or condition of the person using the appliance. Stop using the appliance immediately, and keep the skin cool and moist.
The pain and redness do not subside after use but is this a problem?	▶ If the symptoms subside within 2 days, use the appliance at a level that the symptoms do not appear. If the symptoms persist after 2 days, consult a dermatologist. (See page 12.)
Dirt on the guard cannot be removed.	▶ If the dirt cannot be removed even after cleaning, the light will become weaker and the guard will need to be replaced. Contact the store where you purchased the unit or a service center authorized by Panasonic for repair. (Repair part: WESWH81P3107)

Question	Answer
Does the cartridge confirmation lamp flash when the level selector switch is pressed?	▶ The cartridge confirmation lamp will flash when the lamp cartridge is attached/unattached but this is not a malfunction. Reset the appliance if the lamp cartridge has been replaced. (See page 16.)

Troubleshooting

Problem	Action
	▶ Check that the level indication lamp (blue) has changed from blinking to glowing.
	▶ Press the flash switch perpendicularly against the skin until it flashes. (See page 13.)
The flash window does not flash.	▶ If 10 continuous flashes have been emitted, temporarily move the appliance away from the skin.
	▶ If the cartridge confirmation lamp (red) is blinking, the lamp cartridge replacement period is up. Replace with a new lamp cartridge. (See pages 16 and 17.)
The flash preparation time is long.	▶ The flash preparation time may become longer to prevent the main body from overheating.

Problem	Action
The power turned off after all lamps blinked and the appliance beeped twice.	► Recommended ambient temperature for operation is 5 – 35 °C. The appliance might not be usable outside of the recommended temperature. (See page 15.)
The power turned off after the cartridge confirmation lamp (red) blinked quickly (twice/second) and the appliance emitted a series of beeps.	► The lamp cartridge is not attached correctly. Attach the lamp cartridge correctly. (See page 16.)
The power turned off after the cartridge confirmation lamp (red) blinked.	► The lamp cartridge replacement period is up. Replace with a new lamp cartridge. (See pages 16 and 17.)
The cartridge confirmation lamp (red) glows.	► The lamp cartridge replacement period is approaching. Prepare a new lamp cartridge. (See pages 16 and 17.)
The power turned off during use.	► For safety reasons, the power will turn off automatically if the appliance is not used for approximately 1 minute after the power is turned on. Turn the power on again.
The appliance can only be used a few times after charging.	► The battery has reached the end of its life. (See this page.)

Problem	Action
The appliance does not operate even after pressing the power switch.	► Press and hold the power switch (for approximately 3 seconds) when turning the power on. ► Charge the appliance. (See page 11.)

If the problems still cannot be solved, contact the store where you purchased the unit or a service center authorized by Panasonic for repair.

Battery life

Approximately 3 years.

If the operating time is significantly shorter even after a full charge, the battery has reached the end of its life. (The battery life may change due to usage and storage conditions.) The battery in this appliance is not intended to be replaced by consumers. Have the battery replaced by an authorized service center.

Removing the built-in rechargeable batteries


Remove the built-in rechargeable batteries before disposing of the appliance.

Please make sure that the batteries are disposed of at an officially designated location if there is one.

This figure must only be used when disposing of the appliance, and must not be used to repair it. If you dismantle the appliance yourself, it may cause malfunction.

- Remove the appliance from the AC adaptor.
- If the appliance can be turned on and emit light, use it until the battery is depleted.
- Wait for at least 10 minutes after turning the power off before disassembling the appliance.

- Perform steps ① to ⑥ and lift the battery, and then remove it.
- Take care not to short circuit the positive and negative terminals of the removed battery, and insulate the terminals by applying tape to them.

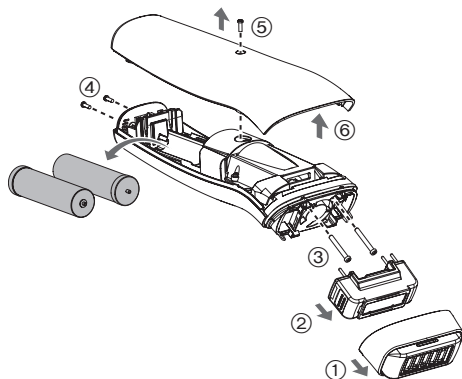For environmental protection and recycling of materials

This appliance contains Li-ion batteries.

Please make sure that the batteries are disposed of at an officially assigned location, if there is one in your country.

Specifications

Power supply	See the name plate on the AC adaptor. (Automatic voltage conversion)
Charging time	Approx. 3 hours

This product is intended for household use only.

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© Panasonic Corporation 2015
Panasonic Corporation